

ENHANCED STRENGTH & POROSITY

Straight & Curved Interbody Fusion Devices



3D Ti Structure with Interconnected Lattice that reduces risk of migration

- Pore Size Range: 100 - 800µm
- Average Pore Size: 600µm

Open Windows, Coupled with Reduced Density (~60% porosity), enhances intraoperative imaging and fusion assessment

Lattice Structure decreases cage stiffness, reducing the risk of subsidence and minimizing stress-shielding





Internal Struts increase overall support and allows for insert and rotate cage positioning

Large Window provides maximum bone graft material

Position Marker ensures correct anatomical placement

Bulleted Leading Edge allows for ease of insertion

Multiple Lordotic Options accommodate varying patient anatomies

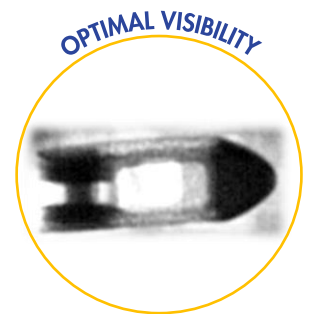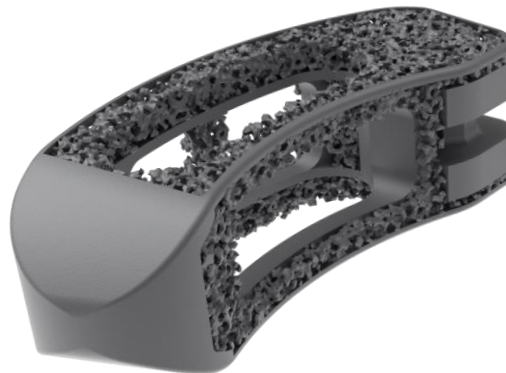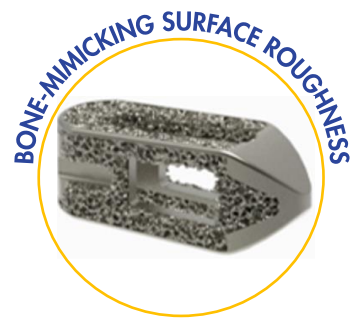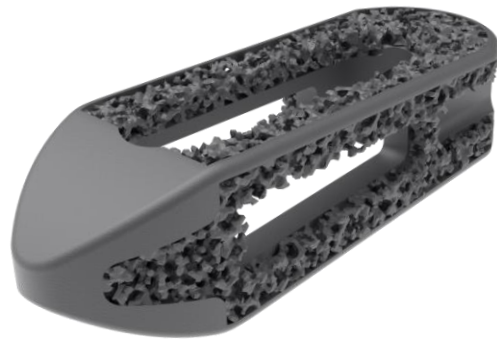

ShurFit® 3D 72-BK-0041 STRAIGHT CAGE KIT



6° & 12° STRAIGHT CAGES (22W x L x H)		6° & 12° STRAIGHT CAGES (24W x L x H)	
P/N	Description (Kit contains 4 of each)	P/N	Description (Kit contains 4 of each)
72-STB-2208	3D Ti Straight Cage; 22 x 08mm 6°	72-STB-2408	3D Ti Straight Cage; 24 x 08mm 6°
72-STB-2209	3D Ti Straight Cage; 22 x 09mm 6°	72-STB-2409	3D Ti Straight Cage; 24 x 09mm 6°
72-STB-2210	3D Ti Straight Cage; 22 x 10mm 6°	72-STB-2410	3D Ti Straight Cage; 24 x 10mm 6°
72-STB-2211	3D Ti Straight Cage; 22 x 11mm 6°	72-STB-2411	3D Ti Straight Cage; 24 x 11mm 6°
72-STB-2212	3D Ti Straight Cage; 22 x 12mm 6°	72-STB-2412	3D Ti Straight Cage; 24 x 12mm 6°
72-STB-2213	3D Ti Straight Cage; 22 x 13mm 6°	72-STB-2413	3D Ti Straight Cage; 24 x 13mm 6°
72-STB-2214	3D Ti Straight Cage; 22 x 14mm 6°	72-STB-2414	3D Ti Straight Cage; 24 x 14mm 6°
72-STC-2209	3D Ti Straight Cage; 22 x 09mm 12°	72-STC-2410	3D Ti Straight Cage; 24 x 10mm 12°
72-STC-2210	3D Ti Straight Cage; 22 x 10mm 12°	72-STC-2411	3D Ti Straight Cage; 24 x 11mm 12°
72-STC-2211	3D Ti Straight Cage; 22 x 11mm 12°	72-STC-2412	3D Ti Straight Cage; 24 x 12mm 12°
72-STC-2212	3D Ti Straight Cage; 22 x 12mm 12°	72-STC-2413	3D Ti Straight Cage; 24 x 13mm 12°
72-STC-2213	3D Ti Straight Cage; 22 x 13mm 12°	72-STC-2414	3D Ti Straight Cage; 24 x 14mm 12°
72-STC-2214	3D Ti Straight Cage; 22 x 14mm 12°		

ShurFit® 3D 72-BK-0040 STRAIGHT CAGE KIT

6° & 12° STRAIGHT CAGES (24W x L x H)		6° & 12° STRAIGHT CAGES (26W x L x H)	
P/N	Description (Kit contains 4 of each)	P/N	Description (Kit contains 4 of each)
72-STB-2408	3D Ti Straight Cage; 24 x 08mm 6°	72-STB-2608	3D Ti Straight Cage; 26 x 08mm 6°
72-STB-2409	3D Ti Straight Cage; 24 x 09mm 6°	72-STB-2609	3D Ti Straight Cage; 26 x 09mm 6°
72-STB-2410	3D Ti Straight Cage; 24 x 10mm 6°	72-STB-2610	3D Ti Straight Cage; 26 x 10mm 6°
72-STB-2411	3D Ti Straight Cage; 24 x 11mm 6°	72-STB-2611	3D Ti Straight Cage; 26 x 11mm 6°
72-STB-2412	3D Ti Straight Cage; 24 x 12mm 6°	72-STB-2612	3D Ti Straight Cage; 26 x 12mm 6°
72-STB-2413	3D Ti Straight Cage; 24 x 13mm 6°	72-STB-2613	3D Ti Straight Cage; 26 x 13mm 6°
72-STB-2414	3D Ti Straight Cage; 24 x 14mm 6°	72-STB-2614	3D Ti Straight Cage; 26 x 14mm 6°
72-STC-2410	3D Ti Straight Cage; 24 x 10mm 12°	72-STC-2610	3D Ti Straight Cage; 26 x 10mm 12°
72-STC-2411	3D Ti Straight Cage; 24 x 11mm 12°	72-STC-2611	3D Ti Straight Cage; 26 x 11mm 12°
72-STC-2412	3D Ti Straight Cage; 24 x 12mm 12°	72-STC-2612	3D Ti Straight Cage; 26 x 12mm 12°
72-STC-2413	3D Ti Straight Cage; 24 x 13mm 12°	72-STC-2613	3D Ti Straight Cage; 26 x 13mm 12°
72-STC-2414	3D Ti Straight Cage; 24 x 14mm 12°	72-STC-2614	3D Ti Straight Cage; 26 x 14mm 12°

ShurFit® 3D 72-BK-0042 STRAIGHT CAGE KIT

6° & 12° STRAIGHT CAGES (26W x L x H)		6° & 12° STRAIGHT CAGES (28W x L x H)	
P/N	Description (Kit contains 4 of each)	P/N	Description (Kit contains 3 of each)
72-STB-2608	3D Ti Straight Cage; 26 x 08mm 6°	72-STB-2808	3D Ti Straight Cage; 28 x 08mm 6°
72-STB-2609	3D Ti Straight Cage; 26 x 09mm 6°	72-STB-2809	3D Ti Straight Cage; 28 x 09mm 6°
72-STB-2610	3D Ti Straight Cage; 26 x 10mm 6°	72-STB-2810	3D Ti Straight Cage; 28 x 10mm 6°
72-STB-2611	3D Ti Straight Cage; 26 x 11mm 6°	72-STB-2811	3D Ti Straight Cage; 28 x 11mm 6°
72-STB-2612	3D Ti Straight Cage; 26 x 12mm 6°	72-STB-2812	3D Ti Straight Cage; 28 x 12mm 6°
72-STB-2613	3D Ti Straight Cage; 26 x 13mm 6°	72-STB-2813	3D Ti Straight Cage; 28 x 13mm 6°
72-STB-2614	3D Ti Straight Cage; 26 x 14mm 6°	72-STB-2814	3D Ti Straight Cage; 28 x 14mm 6°
72-STC-2610	3D Ti Straight Cage; 26 x 10mm 12°	72-STC-2810	3D Ti Straight Cage; 28 x 10mm 12°
72-STC-2611	3D Ti Straight Cage; 26 x 11mm 12°	72-STC-2811	3D Ti Straight Cage; 28 x 11mm 12°
72-STC-2612	3D Ti Straight Cage; 26 x 12mm 12°	72-STC-2812	3D Ti Straight Cage; 28 x 12mm 12°
72-STC-2613	3D Ti Straight Cage; 26 x 13mm 12°	72-STC-2813	3D Ti Straight Cage; 28 x 13mm 12°
72-STC-2614	3D Ti Straight Cage; 26 x 14mm 12°	72-STC-2814	3D Ti Straight Cage; 28 x 14mm 12°

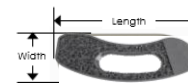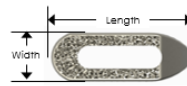

ShurFit® 3D 72-BK-0021 CURVED CAGE KIT

6° & 12° CURVED CAGES (24W x L x H)		6° & 12° CURVED CAGES (28W x L x H)	
P/N	Description (Kit contains 3 of each)	P/N	Description (Kit contains 3 of each)
72-CTA-2406*	3D Ti Curved Cage; 24 x 06mm 0°	72-CTA-2806*	3D Ti Curved Cage; 28 x 06mm 0°
72-CTB-2407	3D Ti Curved Cage; 24 x 07mm 6°	72-CTB-2807	3D Ti Curved Cage; 28 x 07mm 6°
72-CTB-2408	3D Ti Curved Cage; 24 x 08mm 6°	72-CTB-2808	3D Ti Curved Cage; 28 x 08mm 6°
72-CTB-2409	3D Ti Curved Cage; 24 x 09mm 6°	72-CTB-2809	3D Ti Curved Cage; 28 x 09mm 6°
72-CTB-2410	3D Ti Curved Cage; 24 x 10mm 6°	72-CTB-2810	3D Ti Curved Cage; 28 x 10mm 6°
72-CTB-2411	3D Ti Curved Cage; 24 x 11mm 6°	72-CTB-2811	3D Ti Curved Cage; 28 x 11mm 6°
72-CTB-2412	3D Ti Curved Cage; 24 x 12mm 6°	72-CTB-2812	3D Ti Curved Cage; 28 x 12mm 6°
72-CTB-2413	3D Ti Curved Cage; 24 x 13mm 6°	72-CTB-2813	3D Ti Curved Cage; 28 x 13mm 6°
72-CTB-2414	3D Ti Curved Cage; 24 x 14mm 6°	72-CTB-2814	3D Ti Curved Cage; 28 x 14mm 6°
72-CTC-2408	3D Ti Curved Cage; 24 x 08mm 12°	72-CTC-2808	3D Ti Curved Cage; 28 x 08mm 12°
72-CTC-2409	3D Ti Curved Cage; 24 x 09mm 12°	72-CTC-2809	3D Ti Curved Cage; 28 x 09mm 12°
72-CTC-2410	3D Ti Curved Cage; 24 x 10mm 12°	72-CTC-2810	3D Ti Curved Cage; 28 x 10mm 12°
72-CTC-2411	3D Ti Curved Cage; 24 x 11mm 12°	72-CTC-2811	3D Ti Curved Cage; 28 x 11mm 12°
72-CTC-2412	3D Ti Curved Cage; 24 x 12mm 12°	72-CTC-2812	3D Ti Curved Cage; 28 x 12mm 12°
72-CTC-2413	3D Ti Curved Cage; 24 x 13mm 12°	72-CTC-2813	3D Ti Curved Cage; 28 x 13mm 12°
72-CTC-2414	3D Ti Curved Cage; 24 x 14mm 12°	72-CTC-2814	3D Ti Curved Cage; 28 x 14mm 12°

ShurFit® 3D 72-BK-0020 CURVED CAGE KIT

6° & 12° CURVED CAGES (28W x L x H)		6° & 12° CURVED CAGES (32W x L x H)	
P/N	Description (Kit contains 3 of each)	P/N	Description (Kit contains 3 of each)
72-CTA-2806*	3D Ti Curved Cage; 28 x 06mm 0°	72-CTA-3206*	3D Ti Curved Cage; 32 x 06mm 0°
72-CTB-2807	3D Ti Curved Cage; 28 x 07mm 6°	72-CTB-3207	3D Ti Curved Cage; 32 x 07mm 6°
72-CTB-2808	3D Ti Curved Cage; 28 x 08mm 6°	72-CTB-3208	3D Ti Curved Cage; 32 x 08mm 6°
72-CTB-2809	3D Ti Curved Cage; 28 x 09mm 6°	72-CTB-3209	3D Ti Curved Cage; 32 x 09mm 6°
72-CTB-2810	3D Ti Curved Cage; 28 x 10mm 6°	72-CTB-3210	3D Ti Curved Cage; 32 x 10mm 6°
72-CTB-2811	3D Ti Curved Cage; 28 x 11mm 6°	72-CTB-3211	3D Ti Curved Cage; 32 x 11mm 6°
72-CTB-2812	3D Ti Curved Cage; 28 x 12mm 6°	72-CTB-3212	3D Ti Curved Cage; 32 x 12mm 6°
72-CTB-2813	3D Ti Curved Cage; 28 x 13mm 6°	72-CTB-3213	3D Ti Curved Cage; 32 x 13mm 6°
72-CTB-2814	3D Ti Curved Cage; 28 x 14mm 6°	72-CTB-3214	3D Ti Curved Cage; 32 x 14mm 6°
72-CTC-2808	3D Ti Curved Cage; 28 x 08mm 12°	72-CTC-3208	3D Ti Curved Cage; 32 x 08mm 12°
72-CTC-2809	3D Ti Curved Cage; 28 x 09mm 12°	72-CTC-3209	3D Ti Curved Cage; 32 x 09mm 12°
72-CTC-2810	3D Ti Curved Cage; 28 x 10mm 12°	72-CTC-3210	3D Ti Curved Cage; 32 x 10mm 12°
72-CTC-2811	3D Ti Curved Cage; 28 x 11mm 12°	72-CTC-3211	3D Ti Curved Cage; 32 x 11mm 12°
72-CTC-2812	3D Ti Curved Cage; 28 x 12mm 12°	72-CTC-3212	3D Ti Curved Cage; 32 x 12mm 12°
72-CTC-2813	3D Ti Curved Cage; 28 x 13mm 12°	72-CTC-3213	3D Ti Curved Cage; 32 x 13mm 12°
72-CTC-2814	3D Ti Curved Cage; 28 x 14mm 12°	72-CTC-3214	3D Ti Curved Cage; 32 x 14mm 12°

PRECISION SPINE®

www.precisionspineinc.com



*By Request