


PRECISION SPINE
REFORM[®]
HA COATED PEDICLE SCREW

IMMOBILIZATION & STABILIZATION

The Reform[®] HA Coated Pedicle Screw System is designed to meet the varying requirements of degenerative, trauma and deformity procedures. Reform features a cobalt chrome tulip, a titanium triple lead proximally tapered thread to deliver strength, stability and efficiency to all thoracolumbar constructs. The screws have hydroxyapatite coating applied to the threads, which largely consists of calcium and phosphorus, and is designed to enhance the amount of fixation possible between the pedicle screw and the surrounding bone in osteoporotic patients. The screws are supplied sterile, and are made from medical grade titanium or titanium alloy.

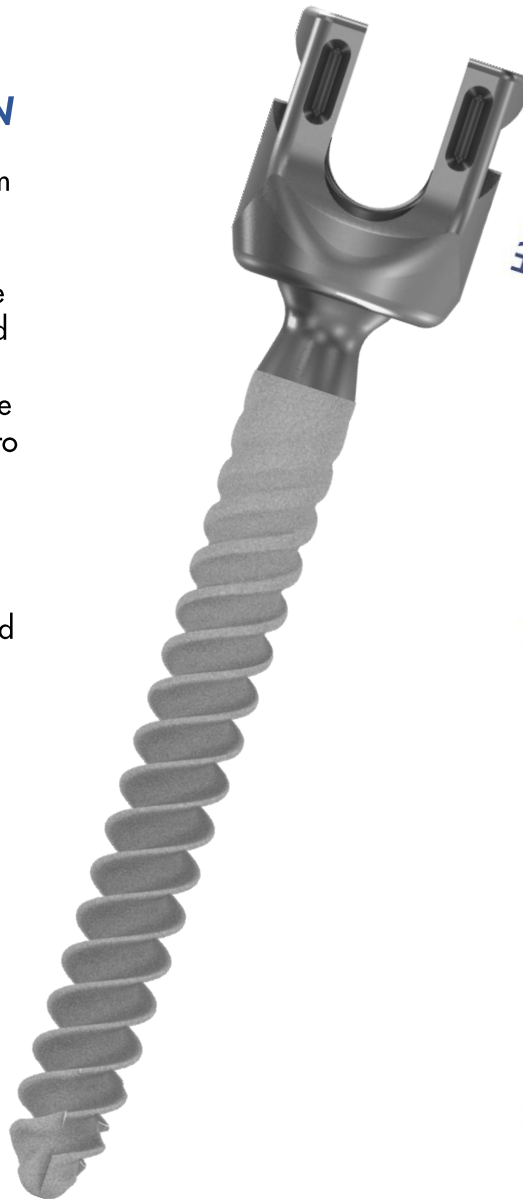


Discover the Difference



REFORM® HA COATED PEDICLE SCREW

Hydroxyapatite Coating largely consists of Calcium and Phosphorus, promoting osteointegration between bone and implant

HA Coating adds approximately .18mm to major diameter and .25mm to minor diameter compared to standard Reform Screws

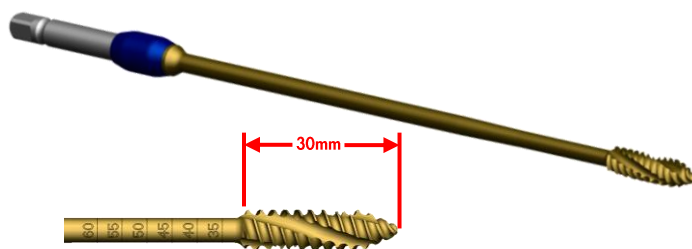

Titanium Nitrate Coated HA Taps are used when high bone density exists

- Standard Reform Taps for Uncoated Screws
 - 0.5mm undersized on the Major Ø
 - Line-to-Line on the Minor Ø



- **Reform HA Taps** for Coated Screws

- 0.5mm undersized on the Major Ø
- Line-to-Line on the Minor Ø



Polyaxial Screws provide 30° angulation in all planes for intraoperative flexibility



Proximal Tapered Thread design increases Bone Screw interface, enhancing pull-out strength



Triple Lead Thread allows for efficient screw delivery

Cobalt Chrome Tulip reduces risk of head splay and increases overall construct stabilization



Hexalobular Drive Fittings on Bone Screw and Locking Cap reduces incidence of toggle and stripping



Reform HA Screw Sizing

5.5mm Diameter	6.5mm Diameter	7.5mm Diameter	8.5mm* Diameter
30mm	30mm	35mm	40mm
35mm	35mm	40mm	45mm
40mm	40mm	45mm	50mm
45mm	45mm	50mm	55mm
50mm	50mm	55mm	60mm
	55mm		70mm
			80mm

Reform Rod Sizing – 5.5mm Ø

Ti Straight	Ti Curved	CoCr Straight	CoCr Curved
80mm	35-80mm (5mm)	80mm	35-80mm (5mm)
100mm	90-120mm (10mm)	100mm	90-120mm (10mm)
120mm		120mm	
200mm		200mm	
300mm		300mm	
400mm		400mm	
500mm		500mm	

*Special Order

Precision Spine, Inc.

2050 Executive Drive, Pearl, MS 39208
 Customer Service: 1.888.241.4773
 Phone: 601.420.4244
 Toll Free: 877.780.4370
 Fax: 601.420.5501
www.precisionspineinc.com

